

Nasal Mask

NM-002-TM

Frame

Simplified design eliminates the need of forehead support. Light and comfortable.

Exhalation Port

Mutes noise and redirects air flow away from bed partner.

Rotating Elbow


360 + 45 degrees rotation prevents pulling of the mask from movement at night.

Cushion



Patented dynamic matching technology minimizes air leakage. Soft and comfortable.

Clip

Easy to attach and detach.

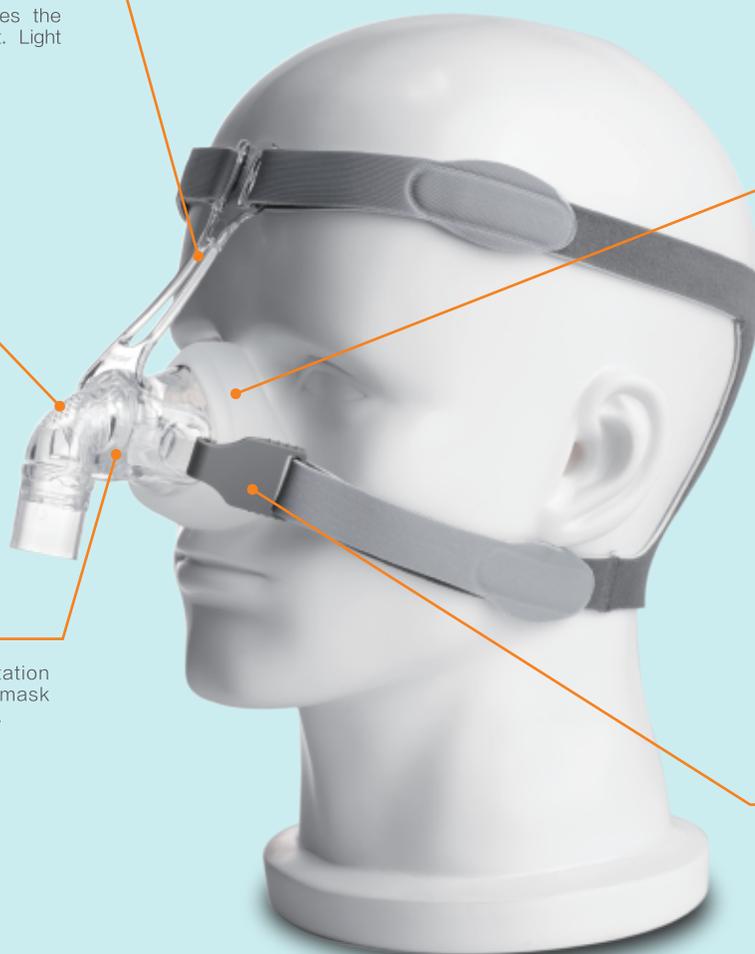

NASAL MASK
PRODUCT STRUCTURE

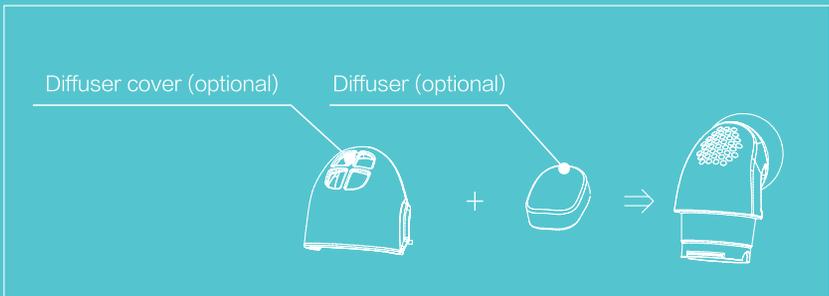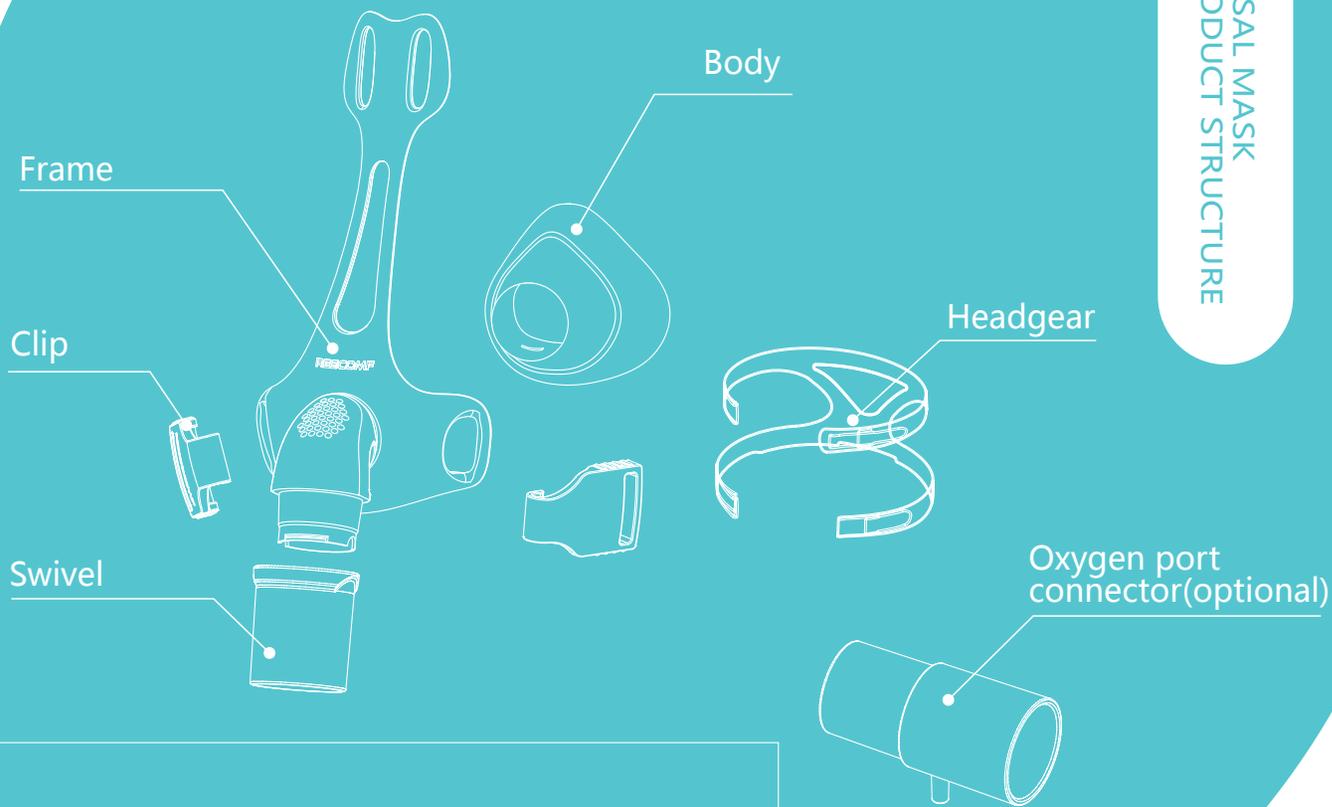

| Comparison with other masks |  |  |  |
|-----------------------------|---|--|---|
| Weight (g) | 44.5 | 88.2 | 89 |
| Size (cm) | 13 × 6.5 × 8.5 | 16 × 7.48 × 8.8 | 16 × 7.8 × 9.0 |
| Noise level (dB) | 20 ± 3 | 22 ± 3 | 22 ± 3 |
| Dead space(cc) | 68 | 126 | 126 |
| Removable cushion | ✓ | × | × |
| Easy cleaning | ✓ | × | × |
| Personalization (option) |  | × | × |

RESCOMF Sales & Marketing Center

Hotline:(86) 400-822-0138

E-mail: info@rescomf.com

Room 2209, Building 5, Huixing Commercial Center, No. 18, Dongsheng Road, Shilong Town, Dongguan City, Guangdong Province, China.

R&D Center : Beijing China